




# ULTIMA A6

## User Manual


### Introduction

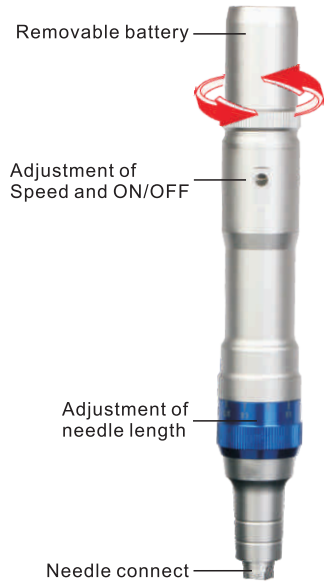

### Parts Name




#### PARTS

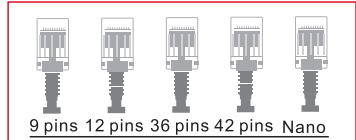

#### MTS

#### Different Cartridges suitable for Device

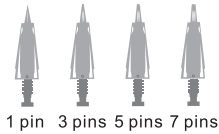

#### PMS

### How to Operate

#### -Battery Mode

##### How to interchange the battery



- (1): unscrew the cover of plug (picture 1)
- (2): replace the battery (picture 2)

##### How to charge the battery



- (1): put the battery into the charger (picture 3)
- (2): Plug in DC adaptor to the jack and then connect power supply (picture 4)



#### Charger

- Indicator light
- The indicator light is red while charging.
- The indicator light is green after fully charged

### How to Operate

#### -Adapter Mode

- (1) Assemble the cover of plug. (picture 5)
- (2) Plug DC adaptor to the jack. (picture 6)
- (3) Plug in DC adaptor to AC outlet.




5

6

#### How to install a cartridges

- (1) Take out needle cartridges in the sterilized packing.
- (2) Fitting round slots 1 into 2 of machine, push needle cartridges horizontally up to the end.
- (3) Lock the cartridge by turning it to the counter clockwise.
- (4) Adjustment of needle length should be done after turn on the machine.
- (5) If you want to remove the needle cartridge, turn it to the clockwise and push out horizontally, fitting round slots (slot1, slot2) together.

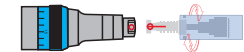

### How to Operate

#### (3) Adjustment of needle length

You can adjust the length of needle by turning the adjustment ring clockwise or counter clockwise. Turning it clockwise, length of it will be longer. Turning it counter clockwise, length of it will be shorter.

#### Needle length adjusting point



Counter Clockwise  
Length of it will be shorter

Clockwise  
Length of it will be longer

### How to Operate

#### 3. Adjustment of speed and needle length

##### (1) Adjustment of Speed and ON/OFF



**Power switch**  
Press 2 seconds to open the device

#### (2) 5 speed level control

Pressing the switch 2 seconds to open the device. It starts from the lowest speed to highest speed when you turn on the power firstly. The speed cycles low to high. Please press the switch 2 seconds when you need to stop the device.

The flicker frequency of light faster, the speed is higher

### Specification

#### 1. Object

Apparatus for applying a Micro-needle Therapy by vibration

#### 2. Performance

A device for puncturing skin-penetrating drugs or cosmetics supports with a micro-needle guided between extended and retracted positions

#### 3. Specification

- 1) Product name: Auto-stamp Motorizedmeso Machine
- 2) Model name: Ultima A6
- 3) Adapter: 5V, 1000mA
- 4) Weight: Main body 62g
- 5) Size: 126\*23(mm)

#### 4. Packing list

Host	1 pcs
Battery	2 pcs
Adapter	1 pcs
Charger	1 pcs
Cartridges	2 pcs 12 pins
Manual	1 pcs

### Maintenance

#### 1. General maintenance

-Keep it with its power off. Wipes its surface out with a dry cloth.  
-Keep it away from directly sunlight or high temperature and humidity.

-Do not leave the measurement, repair or alter of the device to non-professional staff.

#### 2. Routine inspection

-Make sure to inspect the device and adapter before use.

-When you use the device which has not been used for a long time, make sure to check if it works properly with safety.

	Please be sure to read instructions before use and follow as described.
	<b>Do not apply it to the man as follow:</b> 1 Who has acute disease or bloody disease over 38°C 2 Who has allergy and hemophilia 3 Who has melanoma of the face 4 Who has high internal heat 5 Who regularly visits hospital 6 Who takes medicine 7 Who is under 12 years old
	Do not use it on the power other than AC 220V. It may cause fire or malfunction.
	Do not use it near the place with high temperature, dust or humidity (Near the humidifier) it may cause malfunction.
	Do not repair, disassemble or alter the controller and the body. It may cause electric shock.
	Do not touch the device or detach the power cords with wet hands. It may cause electric shock.
	Do not use the stripped cord or propoped cord. It may cause fire (be sure to replace it with another).
	Do not use the propoped cords which has crack or scars. It may cause scars on the skin.
	When you don't use it for a long time, unplug it from the outlet.
	Be sure to have an inspection for test. If any troubles are found, get it repaired.
	No one other than the operator can touch it. It may cause accidents or malfunction.

### Precaution

	Alternating current		Serial No
	Direct current		EU Representative
	Class II equipment		CE mark
	Attention, consult accompanying documents		Manufacturer
	Date of manufacture		Floating applied part

### Warranty

1. We guarantee the products as follow.
2. We offer free repair to our products within 6 months basis (warranty) from the day they are purchased. If any of our products has damages or natural trouble, we offer repairs or replacements.
3. If they are included in the following problems, you have to pay for the repair even during warranty period.
  - Troubles due to a user's careless treatment, voluntary repair or alteration
  - Trouble due to the repairs by a man other than technician from a dealer shop or the customer center appointed by manufacturer
  - Troubles or damages due to shocks while in movement
  - Trouble due to the accessories or options not appointed by manufacturer
  - Trouble due to other natural disaster (fires, salt or water, etc)