




**Dr. pen M8**  
USER MANUAL  
PLEASE READ THE MANUAL BEFORE USE

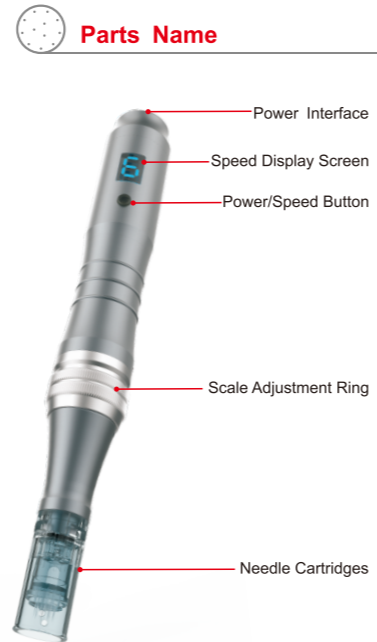


**Parts Name**

**How to operate**

**Wireless Model: Dr pen M8-W**  
The instrument used with both built-in battery or AC power

**Wired Model: Dr pen M8-C**  
The instrument used with with AC power only


- 1.Plug in DC adaptor to the jack
- 2.Plug in DC adaptor to AC outlet
- 3.To charge the battery:Plug the adapter into an outlet and connect it to the power interface

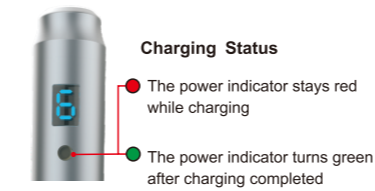

**How to operate**

**Preparation**

- 1.Take out needle cartridge from the sterilized packaging
  - 2.Fitting slot of cartridge into the pen correctly, push needle cartridge horizontally up to the end.
  - 3.Lock the needle cartridge by turning it to the counterclockwise
  - 4.Turn the needle cartridge clockwise and then pull out horizontally if you want to remove it. ( Fig1 )
  - 5.Start the device
- Press the power button for 2-3 seconds to start or shut down the device.

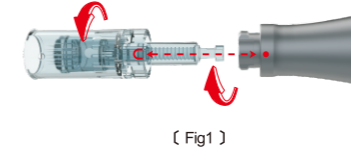

**Attention**  
Before assembling or disassembling the needle cartridge, make sure the U-shaped point of the needle cartridge is aligned with the point of the handpiece

**Attention**  
In order to protect the pen well,do not invert the instrument to prevent the liquid backflow

**How to operate**

**Adjustment of needle length**

Adjust the depth of needle by turning the scale adjustment ring. The depth of needle should be adjusted after start the pen.


Turn it clockwise,needle depth will be shorter  
Turn it counterclockwise,needle depth will be longer



**6 Speeds Adjustment**  
The pen starts from low to high speed  
The speed changes once each time you press the button

**Dr. pen**  
USER MANUAL

**Needle cartridge introduction**

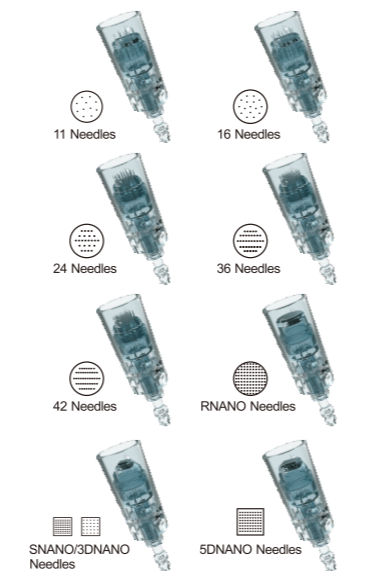

**Specification**

**Product name:**  
**Auto-stamp Motorized Meso Machine**

Model name: Dr pen M8  
Adapter: 5V-1000mA  
Package Weight: about 333g  
Package Size: 169\*123\*56mm

**Packing List**

Item	Quantity
Hand Piece	1 pc
Adapter	1 pc
Manual	1 pc
Power Cable	1 pc
Needle Cartridge	2 pcs

**Precaution**

	Please be sure to read instructions before use and follow as described.
	<b>Do not apply it to the man as follow:</b> 1 Who has acute disease or bloody disease over 38 °C 2 Who has allergy and hemorrhage 3 Who has melasma of the face 4 Who has high internal heat 5 Who regularly visits hospital 6 Who takes medicine 7 Who is under 12 years old
	Do not use it on the power other than AC 220V. It may cause fire or malfunction.
	Do not use it near the place with high temperature, dust or humidity (Near the humidifier).It may cause malfunction.
	Do not repair, disassemble or alter the controller and the body. It may cause malfunction.
	Do not touch the device or detach the power cords with wet hands. It may cause electric shock.
	Do not use the stripped cord or propped cord. It may cause fire(Be sure to replace it with another).
	Do not use the propped cord which has crack or scars. It may cause scars on the skin.
	When you don't use it for a long time, unplug it from the outlet.
	Be sure to have an inspection for test. If any troubles are found, get it repaired.
	No one other than the operator can touch it. It may cause accidents or malfunction.

	Alternating current		Serial No
	Direct current		CE mark
	Attention, consult accompanying documents		Manufacturer
	Date of manufacture		Floating applied part

**Warranty**

The warranty period of device is 6 months. We offer repair to our products from the day they are purchased. If any of our products has damages or natural trouble. We offer repair or replacements

**General maintenance**

Make sure to inspect the device and adapter before use  
Keep it with its power off.Wipe its surface out with a dry cloth

Keep it away from direct sunlight or high temperature and humidity

Do not leave the measurement, repair or alter of the device to non-professional staff

When you use the device which has not been used for a long time, make sure to check if it works properly with safety